Trouble-Shooting:

Problem	Solutions
<ul style="list-style-type: none"> • Not hot enough • No Power • F1 or F2 displayed 	<ul style="list-style-type: none"> • This product may take approximately 15 minutes to acheive selected heat level. • Verify no physical damage to the cord, controller or pad. • Verify tight connection between pad and controller. Plug into a 120V AC outlet. • Verify AC outlet is functional or try an alternate outlet. • If failure continues, discontinue use. See Warranty Information to contact Consumer Services. • After evaluating above and F1 or F2 continues to be dispalyed, seek warranty service.
<ul style="list-style-type: none"> • Control Shuts Off • F1 or F2 displayed 	<ul style="list-style-type: none"> • Auto Off feature is set to turn unit off after selected auto-off time ends. • Re-start unit by pressing On Button. • Verify tight connection between pad and controller.
<ul style="list-style-type: none"> • Pad has a Hot Spot 	<ul style="list-style-type: none"> • Make sure pad is not folded or bunched up in any area.

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5-YEAR LIMITED LIFETIME WARRANTY

Sunbeam Products, Inc. doing business as Jarden Consumer Solutions or if in Canada, Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions (collectively "JCS") warrants that for a period of five years from the date of purchase, this product will be free from defects in material and workmanship. JCS, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty.

This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. JCS dealers, service centers, or retail stores selling JCS products do not have the right to alter, modify or any way change the terms and conditions of this warranty.

This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than JCS or an authorized JCS service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

What are the limits on JCS's Liability?

JCS shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition.

Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty.

JCS disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise.

JCS shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party.

Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How To Obtain Warranty Service**In the U.S.A.**

If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-458-8407 and a convenient service center address will be provided to you.

In Canada

If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-667-8623 and a convenient service center address will be provided to you.

In the U.S.A., this warranty is offered by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions located in Boca Raton, Florida 33431. In Canada, this warranty is offered by Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions, located at 20 B Hereford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.

© 2009 Sunbeam Products, Inc. doing business as Jarden Consumer Solutions. All rights reserved.
Distributed by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions, Boca Raton, Florida 33431.
Printed in China.

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Model 828-511

IMPORTANT NOTES

Please read this entire manual carefully before using this device.
Product may vary from picture.

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www.sunbeam.com

P.N. 137117

IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

DANGER: TO REDUCE RISKS OF BURNS, ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

1. READ ALL INSTRUCTIONS
2. DO NOT USE WHILE SLEEPING.
3. DO NOT USE ON AN INFANT OR ON AN ANIMAL.
4. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, A SLEEPING OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, OR A PERSON WITH POOR BLOOD CIRCULATION.
5. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
6. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING. CHECK SKIN UNDER PAD FREQUENTLY.
7. DO NOT USE IN AN OXYGEN ATMOSPHERE.
8. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
9. DO NOT SIT ON OR AGAINST OR CRUSH PAD-AVOID SHARP FOLDS.
10. NEVER PULL THIS PAD BY THE SUPPLY CORD.
11. DO NOT USE THE CORD AS A HANDLE.
12. THIS UNIT SHOULD NOT BE USED BY CHILDREN WITHOUT ADULT SUPERVISION.
13. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF THE BODY NEEDING HEAT.
14. USE THIS PAD ONLY ON A 110-120 VOLT AC CIRCUIT. UNPLUG WHEN NOT IN USE.
15. DO NOT TAMPER WITH THIS PAD IN ANY WAY.
16. THERE ARE NO USER-SERVICEABLE PARTS. IF FOR ANY REASON THIS PAD DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY FOR CONSUMER SERVICE INFORMATION.
17. DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS, ESPECIALLY ONES THAT CONTAIN HEAT-PRODUCING INGREDIENTS. SKIN BURNS COULD RESULT.
18. DO NOT BEND OR PINCH CORD AND DO NOT USE IF CORD IS DAMAGED.
19. NEVER LEAVE THE APPLIANCE UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
20. LOOP CORD LOOSELY WHEN STORING. TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS.
21. SAVE THESE INSTRUCTIONS.

This appliance has a polarized plug (one blade is wider than the other). This plug fits a polarized outlet only one way; if the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician.
DO NOT ATTEMPT TO DEFEAT THIS SAFETY FEATURE.

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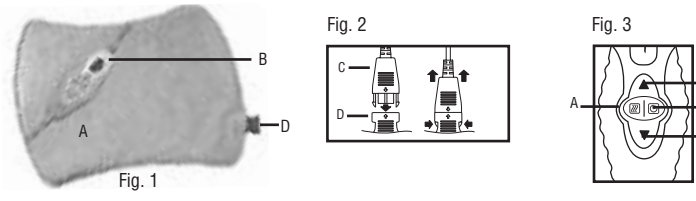
A. Machine washable Flexible Heating Pad
B. LCD Controller
C. Control connector
D. Pad Connector

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HOW TO USE YOUR HEATING PAD

1. Remove all protective wrap or packaging from product before use.
2. Make sure Control Connector is securely connected to the Heating Pad Connector. Push both parts together for a tight fit. Note: The connectors are keyed to fit only one way; align arrow and push together tightly (Figure 2C/D).
3. Plug into 120V AC outlet.
4. To turn on, press On/Off Button (Figure 3A). All control display segments will light up for a few seconds.
5. When turned "On," the LCD controller will automatically default to the Heat mode Selector. "L" will be displayed and indicates that the pad is set to the lowest heat setting.
6. To select a different heat setting from that displayed, press the Up and Down Adjustment Control arrows (Figure 3C). The Heat Mode Selector icon will appear for about 15 seconds. Choose from "L" (Low), 2-7 and "H" (High) settings for your personal comfort.
7. To choose from a range of Auto-off times or a continuous ON setting ("STAY ON"), press the Auto Off Adjustment Selector (Figure 3B). The Auto Off icon will appear for about 15 seconds. Press up and down Adjustment Control Arrows (3C) to select desired Auto off duration. Note: The Auto Off time always defaults to "2 hrs" after the controller is turned off. Auto Off time ranges include: .5 (30 minutes), 1.0 hour, 1.5 hours, 2 hours, 3 hours, 4 hours and "STAY ON."
8. To operate heating pad continuously, without the Auto Off feature, choose the "Stay On" option in step 7. The Stay On icon will remain on display until another Auto-Off adjustment is made or unit is turned off. Please note: BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING. CHECK SKIN UNDER PAD FREQUENTLY. Use extra caution during "Stay On" mode.
9. The controller always defaults back to the Heat Selection mode after either Auto Off or Stay On selection has been made. If an Auto Off time is selected, the controller will turn off automatically after the selected time has elapsed, and the LCD will be turned off after about 15 seconds

USING MOIST HEAT THERAPY

1. Dampen one side of the contoured pad with water by patting the surface with a wet sponge or wet cloth, or use a spray bottle to apply a fine mist.
2. Do not spray or dampen the Pad Connector (Figure D)
DANGER: Do not immerse pad in water or place under running water to dampen surface. Failure to comply may result in personal injury.

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DO NOT DRY CLEAN! DO NOT USE BLEACH!

DO NOT USE WRINGER!

1. Disconnect heating pad from outlet.
2. Detach Pad Connector from Control Connector by pressing tabs on the side of the connector while pulling the two connectors apart.
3. Machine wash pad in cold water on gentle cycle, or hand wash. Tumble dry on low heat or line dry.
4. Do not plug in pad and pad connector until they are completely dry

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